

Topi-CLICK Micro

Revolutionary Way to Accurately & Easily Micro-Dose Topical Medications for Humans & Pets

new

Topi-CLICK[®] MICRO[™]

The Topi-CLICK[®] Micro[™] Difference

Accurate & Consistent Dose[™]

- Topi-CLICK[®] dispensers proven most accurate among others tested by *ARL
- Eliminates user-variability
- Consistent automated dosing throughout treatment

Automated-Metered Dosing

- Audible, Tactile and Visual CLICKs indicate measure of dose at each point of the unique triad shaped base
- One CLICK = 0.05 mL

Targeted Topical Therapy


- Innovative triad shaped dispenser allows for precision handling during application
- Applicator reduces the risk of hand transference to others
- Offers applicator selections for customized dosing



Viscosity Assurance[™]

- Accuracy NOT affected by varying formulations
- Accuracy NOT affected by varying bases

Dispense/CLICK: 0.05 mL Approx. Residual: 0.15 mL
Delivery Volume: 1 mL to 9.75 mL

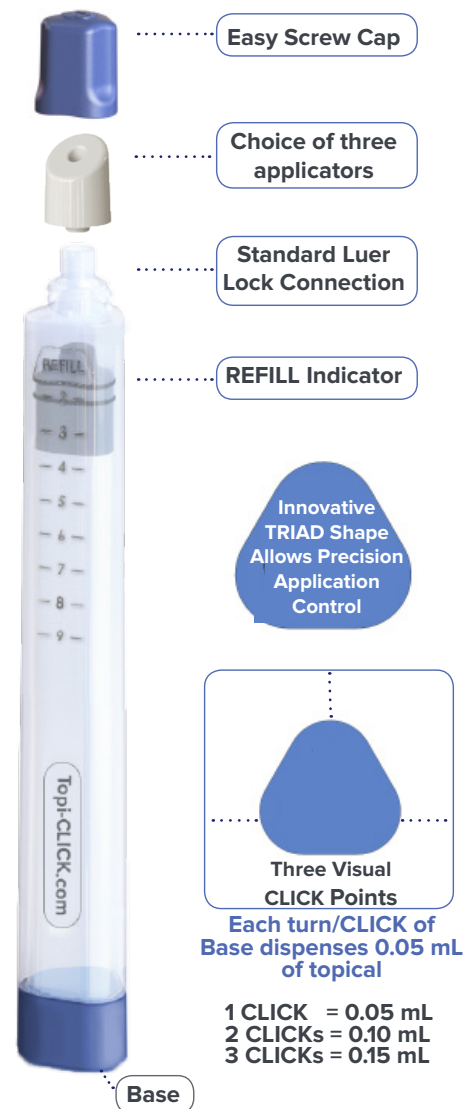

Made of FDA grade materials | Phthalate free | UV block and is Recyclable ♻️.

facebook.com/letcomedical
 linkedin.com/company/letco-medical

460 E. Swedesford Road • Wayne, PA 19808

800.239.5288 | www.letcomed.com

Letco MED
LET US HELP.

Topi-CLICK Micro


User Instructions

- ✓ Automated-metered topical dosing provides an accurate and consistent dose every time
- ✓ Audible, Tactile and Visual CLICKs indicate measure of dose (no need to manually mark dose)
- ✓ Innovative TRIAD shape allows for precision handling
- ✓ Built-in REFILL Badge helps patients reorder medication on-time



STEP 1:

Unscrew the cap of the Topi-CLICK® Micro™ and set aside for recapping after use.

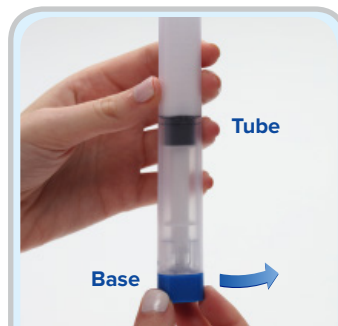

STEP 2:

Holding the Micro™ in an upright position, turn/CLICK the **Base** the number of CLICKs or mLs specified on the prescription label. Use the Audible, Tactile and Visual CLICK points of the triad shaped **Base** as dosing indicators. Be sure to align **Base** with the **Triad Tube** to ensure dose is fully dispensed.



STEP 3:

Apply topical medication with built-in applicator as directed by physician or pharmacist. **Continue rubbing until preparation is completely absorbed into the skin.** **IMPORTANT:** The built-in applicator lessens the risk of transference to unintended areas. No need to use your fingertips to apply medication.



STEP 4:

Screw-on cap after each dose to prevent evaporation. Be sure cap is secured completely.



The Refill Badge™ reminder system helps patients reorder in a timely manner.

When the **TOP** black line of the **Plunger** reaches the embossed line below the word "REFILL" approximately 32 CLICKs are left. Approximately 12 CLICKs remain at the line above the word "REFILL".

Note: Follow instructions on prescription label before use. Topi-CLICK® Micro™ is not re-fillable.

Made of FDA grade materials | Phthalate free | UV block and is Recyclable ♻️.

460 E. Swedesford Road • Wayne, PA 19808

800.239.5288 | www.letcomed.com



Topi-CLICK Micro

Topi-CLICK MICRO™ Filling Instructions Micro Quick-Fill Technology™ Makes Filling Fast & Simple



1. BEFORE FILLING

Prime syringe or EMP Jar to insure only medication (no air) is loaded into the Micro™.

2. FILLING

Note: Graduated Lines are for filling only and not for patients. **Top** black line of the plunger indicates **Fill Volume**. The embossed **Refill Badge** helps patients reorder medication in a timely manner.



Tip-to-tip Filling

Syringe Tip-to-Tip Filling:


Screw Luer syringe to the Luer opening of the Micro™. Requires no extra Luer Lock adaptor.



EMP Jar Quick-Fill Adaptor™

EMP Jar Quick-Fill Adaptor™ Filling:

Screw-on EMP Jar Quick-fill Adaptor™. Fill topical medications directly from mixing jar with this added Micro™ accessory.

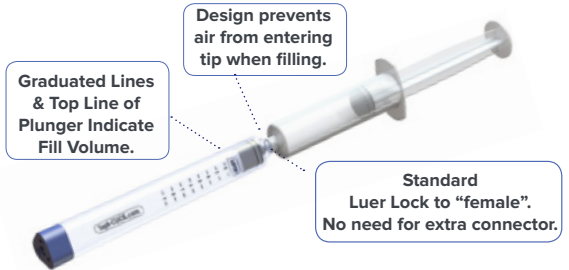

3. APPLICATOR TIP ATTACHMENT



Press on by pushing at a 90 degree angle (important) against a hard surface.

Select one of 3 applicators



Note: Soft Angled Applicator™ placement requires tip to be aligned where lower slant of tip is on the same side as **Refill Badge** of the Micro™ tube.

4. PRIMING MICRO™



Base

Applicator Tip must be attached before priming. Prime the Micro™ by turning/CLICKing the base until contents are dispensed, then wipe off. Be sure to align Micro™ base with the triad tube. Approximate Priming: 0.10 mL or 2 CLICKS.

5. CAP PLACEMENT



Screw-on Cap securely to prevent evaporation.

6. LABEL APPLICATION



Apply label to the bottom side of the Topi-CLICK® Micro™, opposite the "REFILL Badge" indicator.

Dispense/CLICK: 0.05 mL
Delivery Volume: 1 mL to 9.75 mL
Approx. Residual: 0.15 mL

460 E. Swedesford Road • Wayne, PA 19808

800.239.5288 | www.letcomed.com



Topi-CLICK Micro



DoseLogix Laboratory Test Results Prototype 170706-1

Product performed per design specifications. Device is both highly accurate and precise, demonstrating dosage tight variance as well as an average dose being +1.29% greater than the rated. Product performance variance was reduced for two CLICKS versus one CLICK.

Product prototype 170706-1 qualifies inspection and is passed.

Accurate Micro Dosing Designed for Patients and Pets

Topi-CLICK® Micro Lifetime Dosing Accuracy (per 2 CLICKS)	
Average Dose	0.101 mL
Standard Deviation	0.004 mL
Coefficient of Variation	4.46%
Average Dose to Rated Volume <small>(2 CLICKS or 0.10 mL)</small>	+1.29%

Topi-CLICK® Micro Lifetime Dosing Accuracy	
Average Dose	0.051 mL
Standard Deviation	0.004 mL
Coefficient of Variation	8.1%
Average Dose to Rated Volume <small>(0.050 mL)</small>	+1.29%

Dispense/CLICK: 0.05 mL
Dispenser Volume: Up to 9.75 mL
Delivery Volume: Up to 9 mL

Approx. Residual: 0.15 mL
Approx. Priming: 0.10 mL or 2 CLICKS

*Volumes are approximate. As required for compounded prescription pharmaceuticals dispensed in volume dispensers; verify the weight to volume conversions thoroughly.

Made of FDA grade materials | Phthalate free | UV block and is Recyclable 

460 E. Swedesford Road • Wayne, PA 19808

800.239.5288 | www.letcomed.com

